

SIRONA INLAB

HOW TO SEND CONNECT CASES!

October 9, 2014



1. When scanning in the CEREC Connect SW it is the same process as in the CEREC SW. Once you have added the patient's info:
2. Choose whether it is a "Single Restoration" or "Bridge Restoration"
3. Choose the "Tooth #'s"



1. Choose whether it's a:
2. Crown or Inlay/Onlay Full Contour
3. Crown Veneered (PFM, PFG, PFZ)





# 1. Choose what shade you would like



# 1. Choose a stump shade if needed



1. Now that everything has been chosen, you are ready to scan. Remember you can always go back and edit before sending it to the lab.



1. Scan:
2. Lower Jaw
3. Upper Jaw
4. Buccal
5. Add Category if needed:
  - A. (Bio copy, Bio reference, etc.)



# 1. Choose your margins



1. Now you are ready to send to your favorite CEREC Lab:

- A. Username
- B. Password
- C. GO!




# 1. Check your order information



1. Choose Lab
2. Choose return date and time



1. Choose gender
2. Choose age
3. Type additional info
4. Add additional files as needed

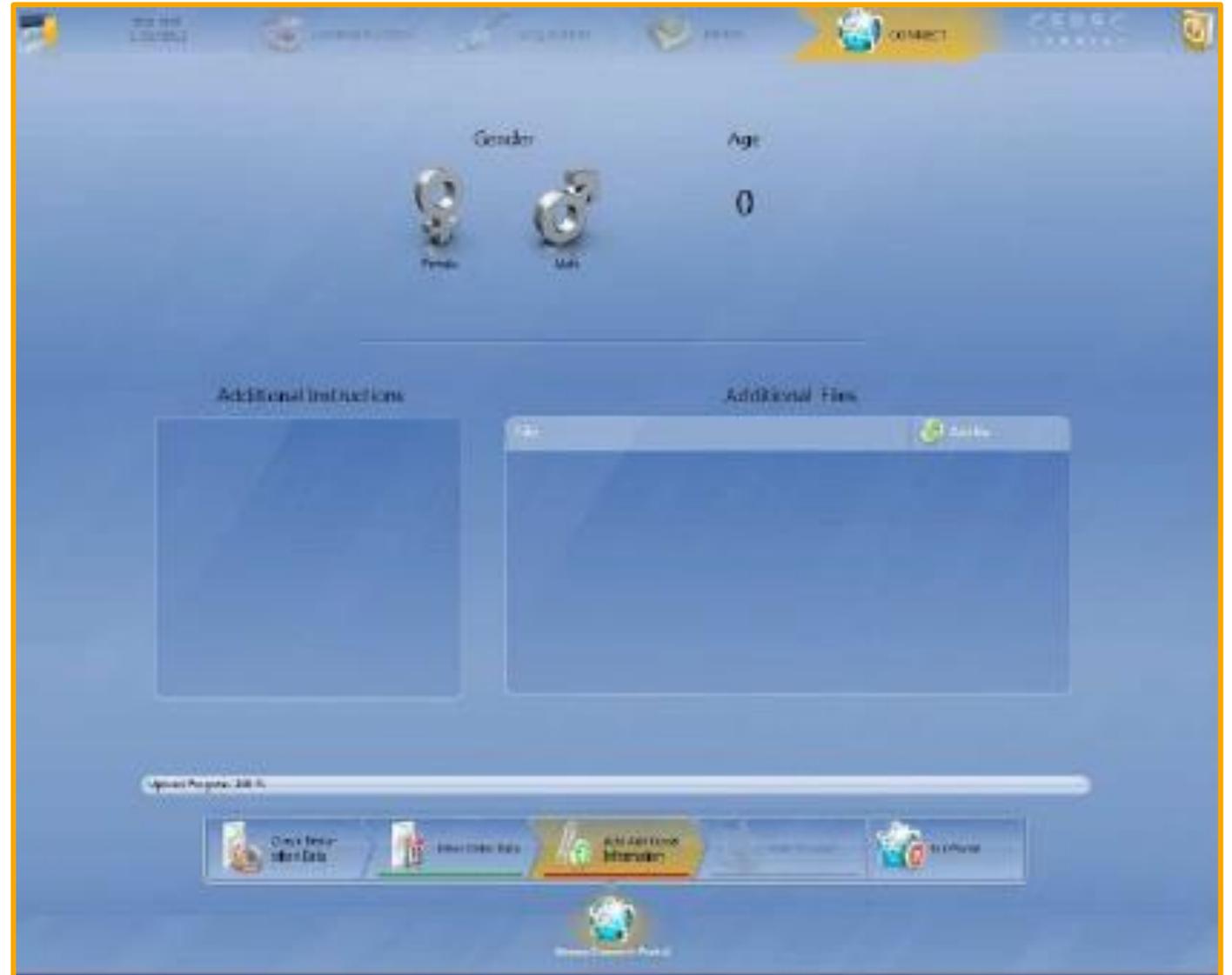

1. At this point the Case has not been sent to the lab yet. It's only in your order cart!
2. Always make sure you "Click the Submit Case" to send to lab! They will never know that a case is pending if this step is forgotten.



1. Once "Submit Case" has been chosen. You must Verify your "Username" and "Password"
2. Click OK. All Done!

The screenshot displays the Sirona website's checkout interface. At the top, the current order number is 100000 (test, test). Below this, there are three main sections: Billing Address, Shipping Address, and Verification. Both the Billing and Shipping addresses are identical, listing Winter Springs Dental Lab at 620 W State Rd 434, Suite 5, Winter Springs, United States. The Verification section contains input fields for Username and Password. Below these sections is the 'Orders in Cart' section, which contains a table with one item: Winter Springs Dental Lab, with a quantity of 100% and a price of \$100.00. At the bottom of the page, there are three buttons: 'View Cart', 'Submit Case', and 'Exit Portal', with 'Submit Case' being the active button.

Order	Part No.	Part	Return Date	Upload Progress		
10000	100,000	Winter Springs Dental Lab	2011/01/01 AM	100% Completed	Cart	Submit